

Switch Adapted Water Gun


USER QUICK GUIDE



Introduction

This is a switch adapted water gun that is compatible with 3.5mm switches. This device may be beneficial to someone with limited hand function that is interested in using a water gun.

Features

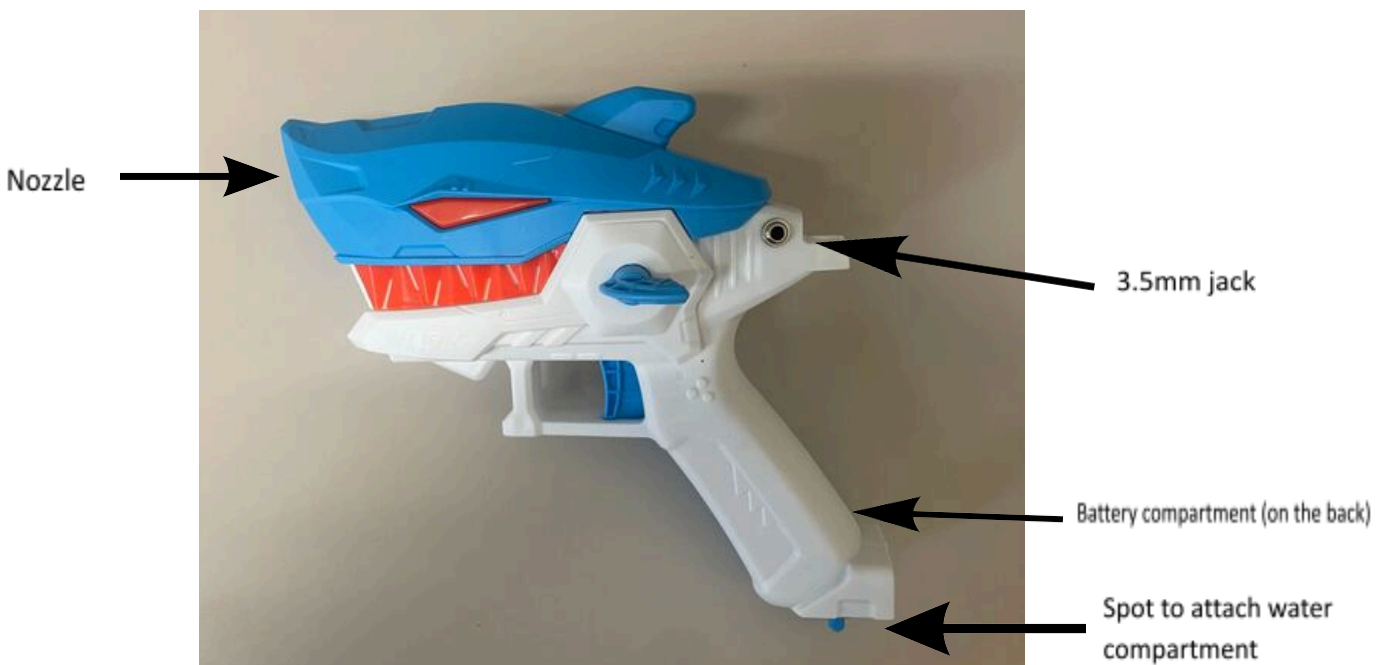

Usage

The Switch Adapted Water Gun can be used by first filling it up with water. Next, plug in a 3.5mm switch into the back of the water gun. Once the switch is pressed, the water gun is activated.

Compatibility

The Switch Adapted Water Gun should be compatible with any 3.5mm switch.

Specifications

Length (in)	Width (in)	Height (in)	Mass (g)
7.95	6.70	7.40	240.4

*The mass given is the mass of the empty water gun.

Cleaning

The device can be wiped down with a damp cloth. When not in use, empty out the water compartment in the gun.



© 2022 by Kerilyn Kennedy – Makers Making Change

This work is licensed under the CC BY SA 4.0 License: <http://creativecommons.org/licenses/by-sa/4.0>

Files available at <https://makersmakingchange.com/project/switch-adapted-water-gun/>

Switch Adapted Water Gun

USER QUICK GUIDE



Install the Battery

The battery compartment, located on the handle of the water gun, can be accessed by removing the screw from it using the provided screwdriver. The battery can be snapped into place with the wires extending from the battery and the water gun. Put the battery in the safety container, place it in the compartment, and screw the compartment shut.

